
# Contour™ Series

Inversion Table Assembly Instructions\*  
For Contour L3 and Contour L5 Models

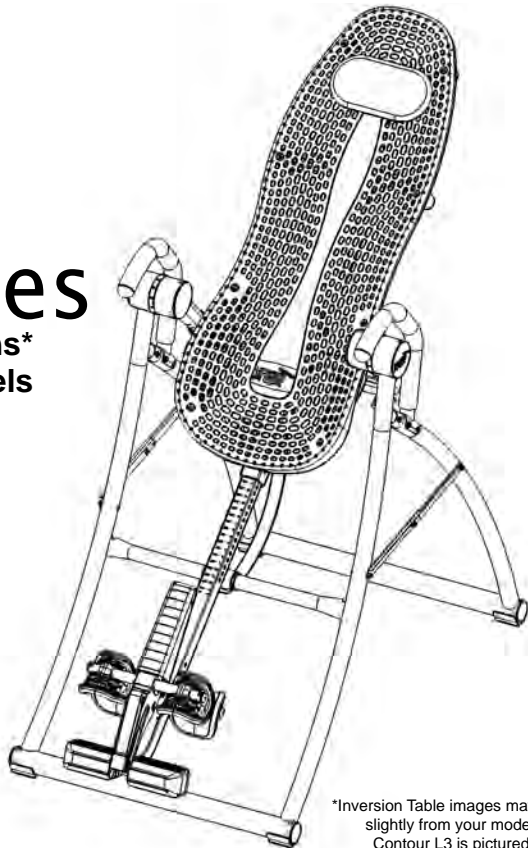\*Inversion Table images may vary slightly from your model. The Contour L3 is pictured here.

**BEFORE YOU BEGIN:** Review all steps before beginning assembly and read all precautions before using the inversion table. Carefully adhere to the Assembly Instructions and Owner's Manual to help ensure safety and product integrity.

**IMPORTANT SAFETY INSTRUCTIONS**  
READ ALL INSTRUCTIONS BEFORE USING THE INVERSION TABLE.

 **WARNING**

**WARNING** - To reduce the risk of injury to persons:

- Read and understand all the instructions, view the Instructional DVD, review all other accompanying documents, and inspect the equipment before using the inversion table. It is your responsibility to familiarize yourself with the proper use of this equipment and the inherent risks of inversion, such as falling on your head or neck, pinching, entrapment, or equipment failure. It is the responsibility of the owner to ensure that all users of the product are fully informed about the proper use of the equipment and all safety precautions.
- Close supervision is necessary when the inversion table is used near children, or by or near invalids or disabled persons.
- Use the inversion table only for its intended use as described in this manual. **DO NOT** use attachments not recommended by the manufacturer.
- **NEVER** drop or insert any object into any opening.
- **DO NOT** use or store product outdoors.
- **DO NOT** use if you are over 6 ft 6 in (198 cm) or over 300 lbs. (136 kg). Structural failure could occur or head/neck may impact the floor during inversion.
- **DO NOT** allow children to use this machine.
- Keep children, bystanders, and pets away from machine while in use.
- Keep body parts, hair, loose clothing and jewelry clear of all moving parts.

**FAILURE TO FOLLOW INSTRUCTIONS AND WARNINGS COULD RESULT IN SERIOUS INJURY OR DEATH.**

**BEFORE YOU BEGIN:** Review all steps before beginning assembly and read all precautions before using the inversion table. Carefully adhere to the Assembly Instructions and Owner's Manual to help ensure safety and product integrity.

## **IMPORTANT SAFETY INSTRUCTIONS**

**READ ALL INSTRUCTIONS BEFORE USING THE INVERSION TABLE.**

### **WARNING**

- The inversion table has no user serviceable parts.
- This product is intended for indoor home use only. **DO NOT** use in any commercial, rental or institutional setting.
- **DO NOT** use the equipment without a licensed physician's approval and a review of the medical contraindications, as noted in the Owner's Manual.
- Failure to assemble and/or use the equipment as directed may void the manufacturer's warranty on this product and could result in injury or death.
- Choose a level surface for assembling and operating the table.
- Follow each step in sequence. **DO NOT** skip ahead.
- Make sure that all fasteners are secure.
- **ALWAYS** test and inspect the table. Make sure the table rotates smoothly to inverted position and back.
- **ALWAYS** replace defective components immediately and/or keep the equipment out of use until repair.

**FAILURE TO FOLLOW INSTRUCTIONS AND WARNINGS COULD RESULT IN SERIOUS INJURY OR DEATH.**

## **SAVE THESE INSTRUCTIONS**


## Which Model Do You Own?

**Important:** Please locate the model of your Teeter inversion table before assembly.

The model can be found on A) the exterior of the box, or B) the serial label located on the back of the table bed.

**Contour L3**

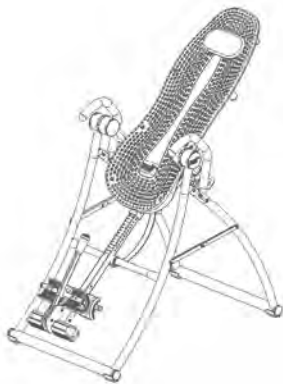**Contour L5**

**My Serial Number:**

\_\_\_\_\_

**Don't forget to register your  
product at [teeter-inversion.com](http://teeter-inversion.com)**

*We hope you enjoy your new Teeter Hang Ups Inversion Table!  
If you have any questions, please contact our Customer Service Department at 1.800.847.0143 or [www.teeter-inversion.com](http://www.teeter-inversion.com).*


# Items for Assembly

Items are not shown to scale. If you have any questions on assembly, contact Customer Service at 1-800-847-0143 or visit our website at [teeter-inversion.com](http://teeter-inversion.com).

ITEMS FOR ASSEMBLY	ITEM NUMBER	ITEMS FOR ASSEMBLY	ITEM NUMBER	ITEMS FOR ASSEMBLY	ITEM NUMBER
<b>A-Frame</b> EZ-Angle™ Tether	NX-1100 F5-1008	<b>Main Shaft</b> with Deluxe EZ-Reach™ Ankle Lock System (Contour L5 model only)	NX-1630-B	<b>Handle Assembly</b> <i>Bolts &amp; Nuts May Come Pre-Assembled</i> Contour Stretch™ Handles (2)	NX-1500 NX-1510 NX-1530
<b>FlexTechnology™ Table Bed</b> Bolts, Washers, Nuts (2) Head Pillow	NX-1300 EP-1128 EP-1105	 with EZ-Reach™ Ankle Lock System (Contour L3 model only)	NX-1630	<b>Tools</b> 10/13 mm Wrenches (2) 5 mm Allen Wrenches (2)	F5-1088 IA-1149
<b>Crossbrace</b> <i>Screw &amp; Post May Come Pre-Assembled</i> Chicago Posts (2) Chicago Screws (2)	NX-1400 F5-0071-A F5-0071-B				

NX-1100

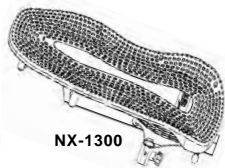NX-1300

EP-1128

EP-1105


NX-1400

F5-0071-A

F5-0071-B

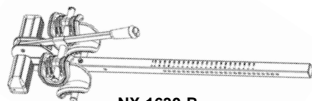NX-1630-B

NX-1630

NX-1500

NX-1510

NX-1530

F5-1088


IA-1149

# Before Beginning

Before reading further, study the drawing below to familiarize yourself with the important components of your new Teeter Hang Ups® inversion table.

1	Head Pillow
2	FlexTechnology™ Table Bed
3	Rotation Adjustment Arm
4	Crossbrace
5	Contour Stretch™ Handles
6	Height-Selector Locking Pin
7	De-Rattler Knob
8	Spreader Arms
9	EZ-Angle™ Tether
10	Crossbar
11	A-Frame
12	Main Shaft
13	Ankle Lock System
14	Ankle Comfort Dial™
15	Stability Feet

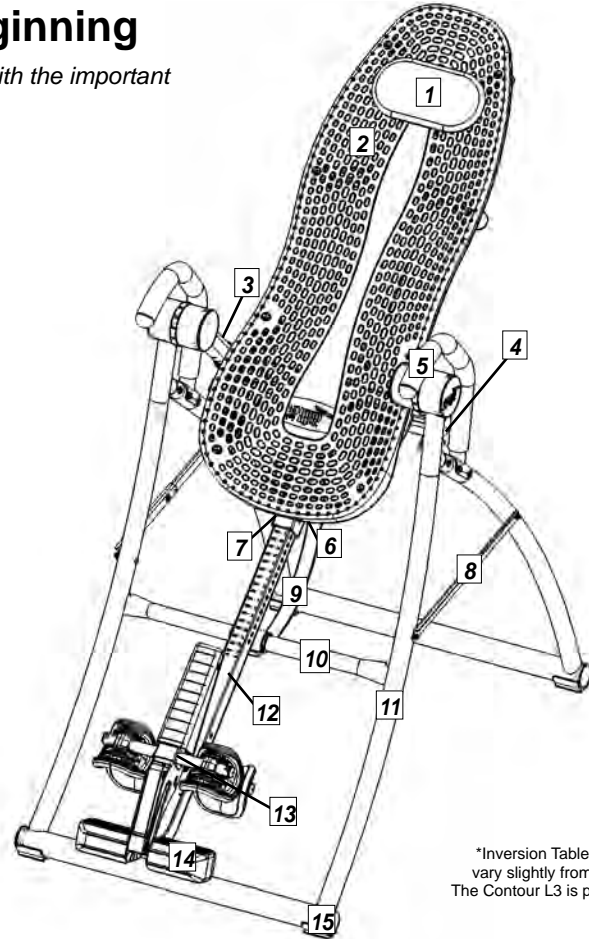\*Inversion Table images may vary slightly from your model. The Contour L3 is pictured here.


# WARNING LABEL PLACEMENT

**Important: Please review all labels and supporting materials before using your inversion table.**

This drawing indicates the locations of the warning labels found on your product.

If a label is missing, illegible or is removed, contact Customer Service at the phone number or website found at the bottom of this page to request a complimentary replacement label.

*NOTE: Image and labels below not shown actual size.*
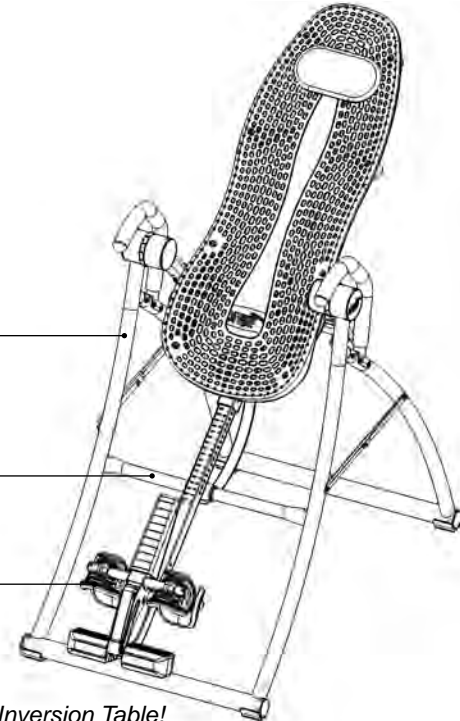
## ▲ WARNING/ AVERTISSEMENT/ ADVERTENCIA

**TIPPING HAZARD:** For upright storage, leave A-Frame open wide enough to remain stable, or secure to the wall to prevent tipping. In households with small children, the table should be stored flat on the floor, not upright.

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## ▲ WARNING

WARNING - To reduce the risk of personal injury or death:

- Read and understand all the instructions before using the inversion table. It is your responsibility to familiarize yourself with the proper use of the equipment and the inherent risks of inversion, such as falling on your head or neck, pinching, entrapment or equipment failure.
- Do not allow children to use the machine.
- Keep children, bystanders, and pets away from the machine while in use.
- Keep body parts, hair, loose clothing and jewelry clear of all moving parts.
- Height/Weight capacity: 4 ft 8 in - 6 ft 6 in (142-198 cm); 300 lbs (136 kgs).
- This product is for consumer, indoor household use only.

Replace Labels and Owner's Manual if Damaged, Illegible, or Removed.  
Teeter, 9902 162nd St. Ct. E., Puyallup, WA 98375  
Toll Free (Phone): 800-847-0143 Web: [www.teeter-inversion.com](http://www.teeter-inversion.com) 48-1228-001-2

## ▲ WARNING AVERTISSEMENT/ ADVERTENCIA

Labels must be properly secured before use.

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
*We hope you enjoy your new Teeter Hang Ups Inversion Table!  
If you have any questions, please contact our Customer Service Department at 1.800.847.0143 or [www.teeter-inversion.com](http://www.teeter-inversion.com).*


# Step 1

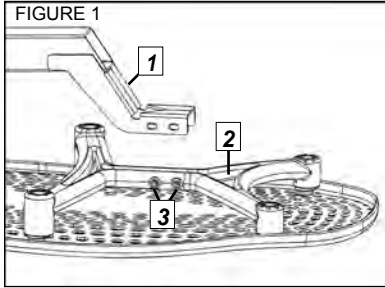- 1 Support Beam
- 2 Upper Support Arms
- 3 Rubber Spacers


## Attach the Support Beam to the Table Bed (NX-1300)

- Place the Table Bed face down on the floor. Push down on the Support Beam so that the two holes align evenly with the holes at the base of the Upper Support Arms (Figure 1).
- Insert the two bolts (EP-1128) into the open holes (Figure 1).
- Fasten each bolt with a washer and nut (Figure 2).
- Use the Allen Wrench (IA-1149) to steady each bolt and tighten the nuts using the 10/13 mm open end Wrench provided (F5-1088) (Figure 2).

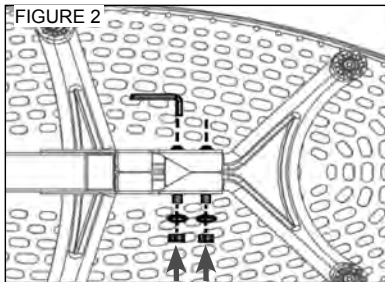
# Step 2


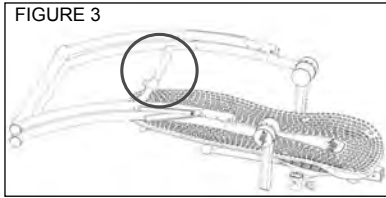
## Assemble the Table Bed (NX-1300) to the A-Frame (NX-1100)

- Position the A-Frame on the floor with the crossbar pointing down, so that the frame is resting on the Rotation Adjustment Arms (Figure 3).
- With the Table Bed face up and the top of the bed pointing to the A-Frame feet, slide it on the floor between legs of the A-Frame.
- Loosen the “de-rattler” function of the Pinlock Knobs by turning each counterclockwise (Figure 4).
- Align and insert the Rotation Adjustment Arms into the Support Housings (Figure 5).
- Pull out both spring-loaded Pinlock Knobs and then pull up to align the Support Housing with one of the marked settings on the Rotation Adjustment Arms (Figure 6).
- Release the knobs so that they engage a hole setting (Figure 6). Should the knobs be locked in the “out” position, simply rotate clockwise to release (Figure 7).

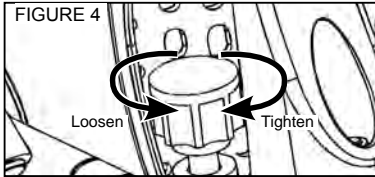
Counterclockwise to Loosen, Clockwise to Tighten

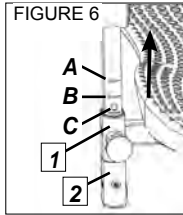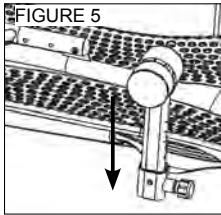- 1 Support Housing
- 2 Rotation Adjustment Arms


*NOTE: Refer to the Owner's Manual for an explanation of the hole settings. If you are unsure, use Setting C to start.*

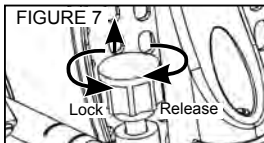Counterclockwise to Lock, Clockwise to Release

## Step 3

FIGURE 8

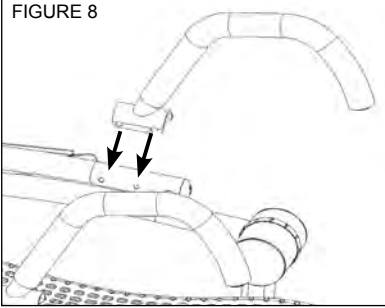
### Attach the Handles (NX-1500) to the A-Frame (NX-1100)

*Bolts and Nuts may come pre-assembled in Handles.*


- To assemble the Handles, leave the inversion table in this position (Step 2). Remove the pre-assembled bolts from the Handles.
- Place each Handle fitting over the A-Frame legs, aligning with the pre-drilled holes (Figure 8).
- Insert the bolts and secure with the nuts, tightening with the Allen Wrenches provided (Figure 9).

FIGURE 9
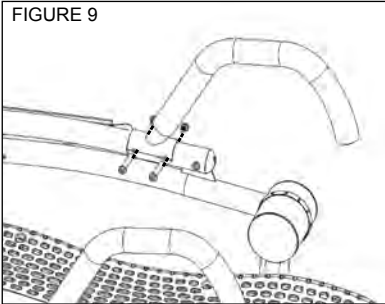
## Step 4

FIGURE 10
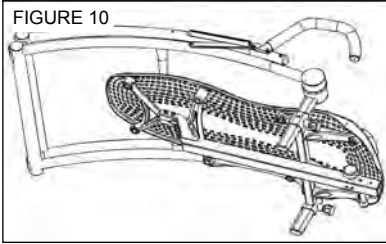
### Stand the Inversion Table Upright

- To stand the inversion table upright, first move the A-Frame onto its side (Figure 10), and then carefully lift it up onto its feet (Figure 11).
- Use the handles to open the A-Frame wide enough so that it is stable and can stand upright on its own (Figure 12).
- Leave the table in this position for Step 5 of the assembly.

FIGURE 11


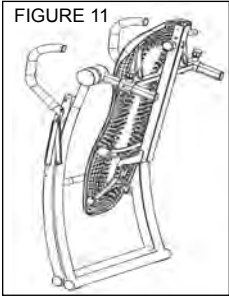
FIGURE 12
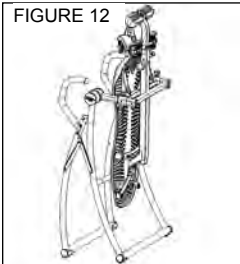
## Step 5

FIGURE 13

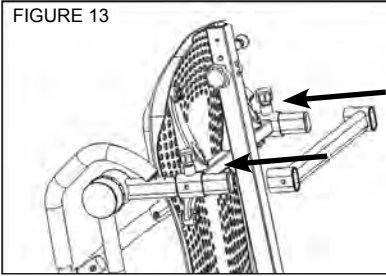
FIGURE 14


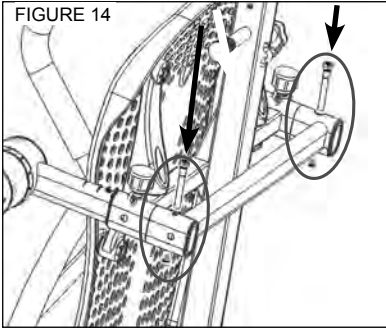
FIGURE 15

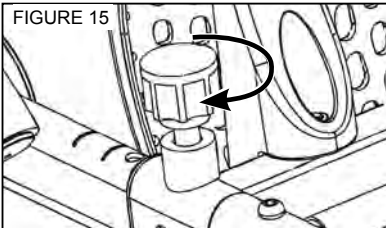Turn Clockwise to Tighten the De-Rattler Function



### Attach the Crossbrace (NX-1400) to the Rotation Adjustment Arms

- Remove the pre-assembled Chicago Screw Posts from the Crossbrace (NX-1400).
- Align the holes in the Crossbrace with the last hole on the Rotation Adjustment Arms (Figure 13).
- Insert the Chicago Posts (F5-0071-A) and secure with the Chicago Screws (F5-0071-B), tightening with the Allen Wrenches (IA-1149) provided (Figure 14).


*NOTE: Due to the tapered style of the Chicago Post, it may not sit flush against the Crossbrace.*

- Turn the Pinlock Knobs clockwise to re-engage the “de-rattler” function (Figure 15).

# Step 6

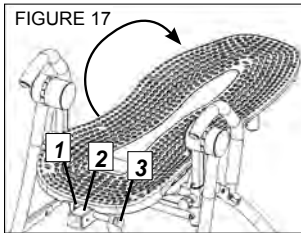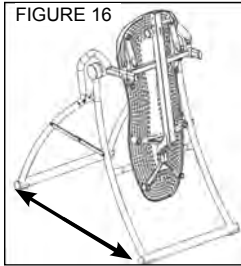- 1 De-Rattler Knob
- 2 Main Shaft Housing
- 3 Height-Selector Locking Pin

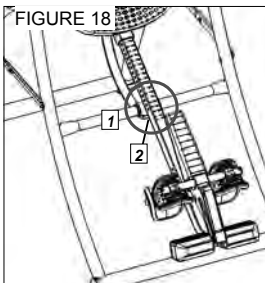- 1 Crossbar
- 2 Proper Resting Position

Depending on the model you have purchased, the items in this step will vary. The model of your product can be found on the Serial Label, located on the back of the Table Bed.

## Insert the Main Shaft into the Table Bed (NX-1300)

**Contour L3:** Use the Main Shaft with the EZ-Reach™ Ratchet Ankle Lock System (NX-1630).

**Contour L5:** Use the Main Shaft with the Deluxe EZ-Reach™ Ankle Lock System (NX-1630-B).

- Ensure that the A-Frame is fully open and locked at the Spreader Arms (Figure 16).
- Rotate the Table Bed into its use position (Figure 17).
- Loosen the Main Shaft De-Rattler Knob at the base of the Main Shaft Housing.

*NOTE: The De-Rattler Knob is an optional feature that prevents a slight shift when you invert. Product security is not compromised if you choose not to tighten the knob.*

- With the Height Adjustment Settings on the Main Shaft facing up, slide the end of the Main Shaft into the Main Shaft Housing at the base of the Table Bed.
- Pull out the Height-Selector Locking Pin to allow the Main Shaft to slide in further and release in the desired height setting. Refer to the Owner's Manual for more information on selecting your height setting.
- The Main Shaft **MUST REST** against the Crossbar of the A-Frame (Figure 18).

*NOTE: The Crossbar prevents the Table Bed from rotating forward when the user steps on the Ankle Comfort Dial. If the Main Shaft does not rest on the Crossbar as shown in Figure 18, then the Table Bed has been assembled backwards onto the A-Frame. This must be corrected before use.*

# Step 7

FIGURE 19


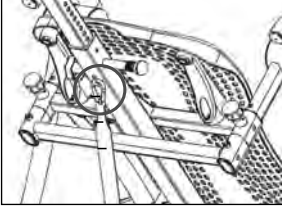
FIGURE 20
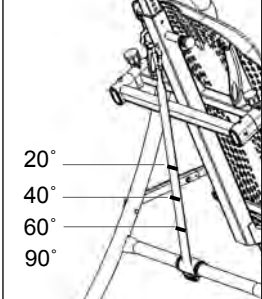
## Attach the EZ-Angle™ Tether Strap (F5-1008) to Limit the Degree of Rotation

- The Tether will come pre-assembled to the A-Frame.
- Unfold the adjustable Tether and clip it to the U-Bar on the underside of the Table Bed (Figure 19).
- The EZ-Angle™ Tether has color coded angles to assist in varied angles of inversion for your experience level (Figure 20):

20° / Twenty Degrees – adjust the strap so the green stripe is in the center of the sliding buckle for inversion to a gentle twenty-degree angle.

40° / Forty Degrees – adjust the strap so the orange stripe is in the center of the sliding buckle for inversion to a moderate forty-degree angle.


60° / Sixty Degrees – adjust the strap so the red stripe is in the center of the sliding buckle for inversion to a more advanced sixty-degree angle. This is the angle when full decompression is realized.

90° / Ninety Degrees – Remove the tether strap clip from the base of the table bed loop to enable rotation to full inversion.

## Step 8


### Attach the Head Pillow (EP-1105) and Owner's Manual

- Attach the Head Pillow by securing the Velcro straps through the holes in the Table Bed – the position of the pillow can be adjusted depending on the user (Figure 21).
- Hang the Owner's Manual over one of the handles for future reference (Figure 21).
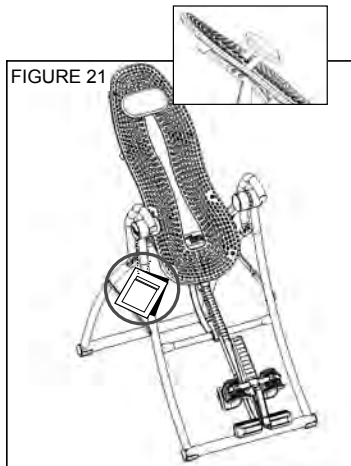
### Before Use

- Test the Inversion Table by hand for smooth rotation (Figure 22).
- Ensure all fasteners are secure.
- Please complete the warranty registration online at [www.teeter-inversion.com](http://www.teeter-inversion.com).
- For your reference, the serial number can be found on the back of the Table Bed.
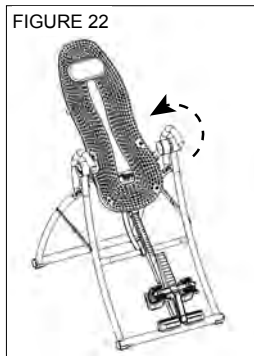
## ⚠ WARNING

Read the Owner's Manual thoroughly before using your Teeter Hang Ups Inversion Table. Improper settings could result in serious injury or death!

# Adjustments

FIGURE 23



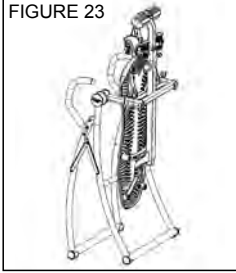
FIGURE 24A

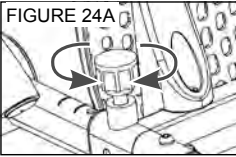Counterclockwise to Loosen,  
Clockwise to Tighten.

FIGURE 24B

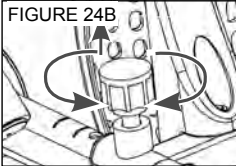Counterclockwise to Lock,  
Clockwise to Release.

**Option:** instead of rotating to lock pins in the “out” position, simply pull the knobs up, slide, and release into your desired setting.


## Changing the Rotation Adjustment Setting

- Loosen the De-Rattler Knob.
  - Shorten the Main Shaft and rotate the Table Bed opposite from use into its storage position. Leave the A-Frame open and locked at the spreader arms (Figure 23).
  - The Pinlock Knobs serve two functions:
    - As a “de-rattler” to prevent excess movement of the support arms (Figure 24A)
    - As a “spring-loaded pinlock” to engage the Rotation Adjustment Setting (Figure 24B).
  - First, loosen the “de-rattler” function of the Pinlock Knobs by rotating them counterclockwise. Then pull the knobs up to disengage the pin and rotate again counterclockwise to lock them in the “out” position.
  - Move the Support Housing along the Rotation Adjustment Arms to the desired setting using the line markings as a guide. Rotate the pin clockwise to release the lock and allow the pin to engage a hole. You will hear a click when the pin is properly engaged.
- NOTE: Engage the Pinlock Knobs in the same setting for each of the Rotation Adjustment Arms.*
- Turn the Pinlock Knobs clockwise again to re-engage the De-Rattler.
  - Rotate the Table Bed back to its use position and extend the Main Shaft to your preferred setting.

# Maintenance / Storage

FIGURE 24
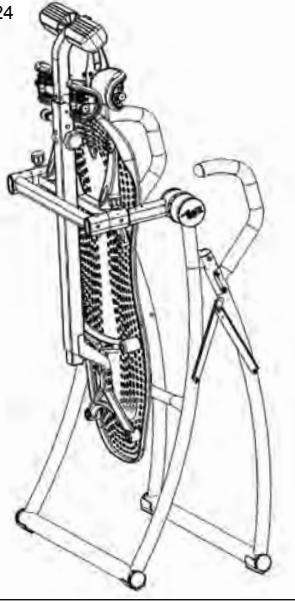
## Maintenance

- To clean the Inversion Table, wipe down with a damp cloth. Do not use abrasive cleaners or solvents.

## Storage

- Remove the Angle Tether from the U-Bar located on the underside of the Table Bed.
- Loosen the Main Shaft De-Rattler Knob.
- Pull the height-selector locking pin, slide the Main Shaft in all the way, and engage the pin in the storage setting.
- Rotate the table opposite from use until the table has turned 180 degrees and rests against the Crossbar on the A-Frame.
- Pull up on the Spreader Arms to fold the A-Frame (Figure 24), leaving the A-Frame legs open to a width of at least 12" for stability.

*NOTE: This operation may pinch fingers if not done slowly and carefully.*

## WARNING

Tipping Hazard: For upright storage, leave A-Frame open wide enough to remain stable, or secure to the wall to prevent tipping. In households with small children, the inversion table should be stored flat on the floor, not upright.

## CONNECT WITH US!

[facebook.com/teeterhangups](https://facebook.com/teeterhangups)

[teeter-inversion.com/blog](http://teeter-inversion.com/blog)

[twitter.com/teeterhangups](https://twitter.com/teeterhangups)

[youtube.com/teetertv](https://youtube.com/teetertv)

If you have any problems assembling the equipment, or questions about its use, please contact customer service at:

USA: Teeter  
Toll Free (Phone) 800-847-0143  
(Fax) 800-847-0188  
[info@teeter-inversion.com](mailto:info@teeter-inversion.com)  
[www.teeter-inversion.com](http://www.teeter-inversion.com)

International: Inversion International, Ltd.  
(Phone) +1-242-362-1001  
(Fax) +1-242-362-1002  
[info@InversionInternational.com](mailto:info@InversionInternational.com)  
[www.InversionInternational.com](http://www.InversionInternational.com)

U.S. and Foreign Patents Pending. Teeter Hang Ups is a registered trademark of Teeter.  
Specifications subject to change without notice.  
© COPYRIGHT 2012 Teeter. International Law Prohibits Any Copying.

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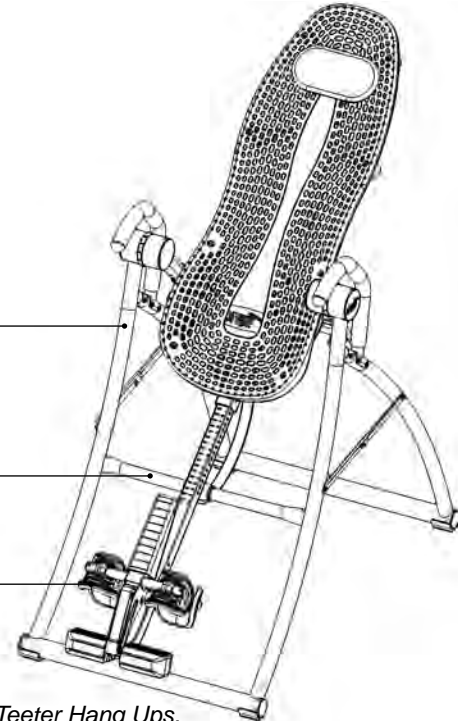
**▲ WARNING/ AVERTISSEMENT/ ADVERTENCIA**

**TIPPING HAZARD:** For upright storage, leave A-Frame open wide enough to remain stable, or secure to the wall to prevent tipping. In households with small children, this table should be anchored to the floor, not upright.

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**▲ WARNING**

**WARNING -** To reduce the risk of personal injury or death:

- Read and understand all the instructions before using the inversion table.
- It is your responsibility to familiarize yourself with the proper use of the equipment and the inherent risks of inversion, such as falling on your head or neck, pinching, entrapment or equipment failure.
- Do not allow children to use the machine.
- Keep children, bystanders, and pets away from the machine while in use.
- Keep body parts, hair, loose clothing and jewelry clear of all moving parts.
- Height/Weight capacity: 4 ft 8 in - 6 ft 6 in (142-198 cm); 300 lbs (136 kgs).
- This product is for consumer, indoor household use only.

**Replace Labels and Owner's Manual if Damaged, Ilegible, or Removed.**

Teeter, 9902 162nd St. Ct. E., Puyallup, WA 98375  
Toll Free (Phone): 800-847-0143 Web: [www.teeter-inversion.com](http://www.teeter-inversion.com) SP-132 9811-3

**▲ WARNING  
AVERTISSEMENT/ ADVERTENCIA**

**ATTENTION:** Read the property manual before use.

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## ⚠ WARNING

**WARNING - To reduce the risk of personal injury or death:**

• Read and understand all the instructions before using the inversion table.  
It is your responsibility to familiarize yourself with the proper use of the equipment and the inherent risks of inversion, such as falling on your head or neck, pinching, entrapment or equipment failure.

- Do not allow children to use the machine.
- Keep children, bystanders, and pets away from the machine while in use.
- Keep body parts, hair, loose clothing and jewelry clear of all moving parts.
- Height/Weight capacity: 4 ft 8 in - 6 ft 6 in (142-198 cm); 300 lbs (136 kgs).
- This product is for consumer, indoor household use only.

**Replace Labels and Owner's Manual if Damaged, Illegible, or Removed.**  
Teeter, 9902 162nd St. Ct. E., Puyallup, WA 98375

Toll Free (Phone): 800-847-0143 Web: [www.teeter-inversion.com](http://www.teeter-inversion.com) 24-107-0013-01

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© COPYRIGHT 2012 Teeter. La Loi Internationale Interdit Toute Copie.